



MEDBOX PLUS PACKAGING FOR TEMPERATURE CRITICAL SHIPMENTS

Specifications

To be used in road and air freight temperature critical shipments, packed in dry ice principally clinical trials (conform ADR and IATA requirements).

	<p>THE BENEFITS OF MEDBOX PLUS:</p> <ul style="list-style-type: none"> ❖ Professional appearance ❖ Complete solution in one pack. ❖ Easy handling. ❖ Possibility of customization. ❖ Regulatory compliance. ❖ Expansion with RFID (track and tracing) is possible.
	<p>PARTS OF A COMPLET MEDBOX PLUS:</p> <ul style="list-style-type: none"> ❖ 1 white cardboard box with inscription. ❖ 1 approved EPS-box with unique fixation ❖ 1 Easygrip bag, containing the parts of the kit ❖ 3 Leak-proof 95kPa Lab-Sealbag® suitable for dry ice ❖ 3 patented cardboard fixing clamps. ❖ 3 absorbent sheets (12x14cm) ❖ 1 instruction manual

DIMENSIONS, FEATURES:


- ❖ Dimensions of the Medbox outer carton (L x W x H): 380 x 380 x 400 mm.
- ❖ Dimensions of the inner EPS-boksz (L x W x H): 300 x 300 x 280 mm.
- ❖ The Medbox can contain up to 16 kg of dry ice.
- ❖ Keeps your shipment frozen up to 100 hours (tested and validated; figure based upon 13.3 kg dry ice).
- ❖ Designed to contain a maximum of 120 pcs 3.5 ml vials.

All statements, technical information and recommendations herein are based on tests the Manufacturer believes to be reliable, but the accuracy and completeness thereof is not guaranteed and the following is made in lieu of all warranties, expressed or implied, including the implied warranties of merchantability and fitness for purpose: Seller's and Manufacturer's only obligation shall be to replace such quantity of the product proved to be defective. Before using, User shall determine the suitability of the product for its intended use, and User assumes all risk and liability whatsoever in connection therewith. Seller reserves the right to improve the product and accordingly change the specification without notice. Neither Seller or Manufacturer shall be liable either in tort or in contract for any loss or damage, direct, incidental or consequential, arising out of the use of or the inability to use the product. No statement or recommendation not contained herein shall have any force or effect unless in an agreement signed by officers of or Seller and Manufacturer.

USER MANUAL




1. Fill the Sealbag® with absorption pouch.. Do not put dry ice in the Sealbag !

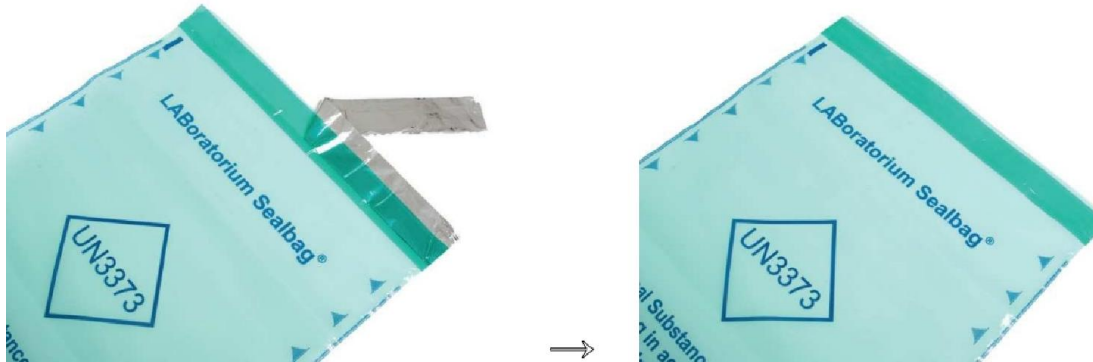

2. Remove the release liner and close the Sealbag®.



3. Fixate the Sealbag® between the cardboard strip and fasten this into the EPS box.

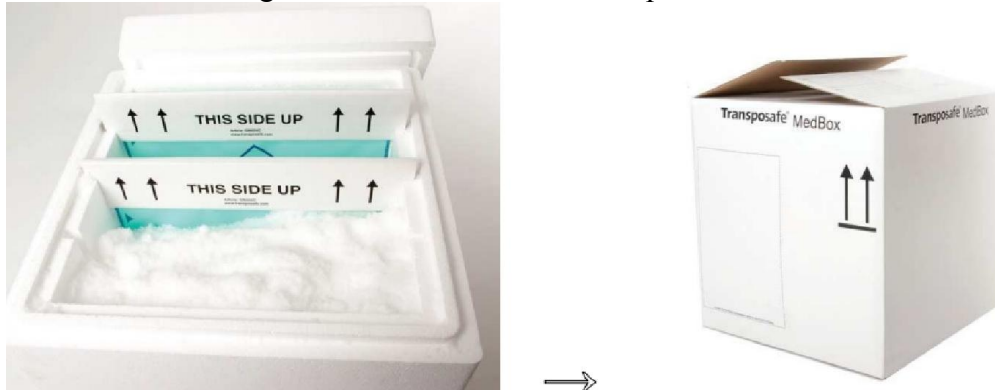

4. Fill the EPS box with dry ice.